

For products: AXM06, AXM09, AXM12, AXM15, AXM18, AXM21, AXM24

Intended Purpose – Portable (folding) ramp intended to assist wheelchair users and persons with reduced mobility to overcome changes in level and environmental barriers.

BEFORE USE


- Retain all packaging
- Read the instructions
- Ensure the surface is stable
- Capacity: check ramp label



Scan for more product information

TO DEPLOY

1. Carry the ramp, using both handles together (see fig 1).
2. Check ramp is long enough to cover the step or rise.
3. The ramps are designed for use between 1:4 and 1:20 gradient (dependent on safe application and use requirements).
4. Place the closed ramp in position with the top lip fully on the upper level and the hinge central to the chosen route (see fig 2). The unique design of the end lips holds the ramp in position.
5. Open the ramp using the 2 handles (see fig 3 & 4). Both hands can be used.
6. Ensure the ramps do not become dislodged in use. Check ramp is firmly located and stable each time it is used.

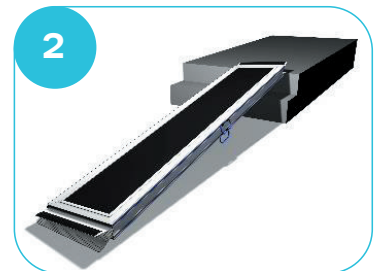

Indication – Impairment of mobility relevant to everyday life caused by steps and / or thresholds

MAINTENANCE

If the ramps become damaged or slippery they should not be used.

WARRANTY

- This product carries a comprehensive 2-year warranty for normal use; conditions apply.
- Use other than published or intended will void the warranty.
- Any modification or repair by anyone other than the manufacturer will void the warranty.
- Removal of the batch or serial number will void the warranty.
- Expected Service Life: 10 years from date of manufacture.
- Any serious incident involving this device should be reported to the manufacturer.

STORAGE

Suitable for indoor and outdoor use but must not be stored outdoors.

RECYCLING

If you no longer require this product, please liaise with a metals recycling contractor. If this product was loaned to you by a local authority, please contact them to discuss its removal.

APPLICATION RISK AND CONTRAINDICATIONS

Do not use if damaged.
Do not exceed the stated maximum load capacity.
Do not use if the ramp cannot be securely positioned.
Do not use on unstable, icy, oily or excessively slippery surfaces.
Ensure the ramp length is appropriate for the height being overcome.
Failure to follow these instructions may result in injury or damage to property.

CLEANING INSTRUCTIONS – The ramp may be cleaned with household cleaning products or disinfectants. Do not use bleach, strong acids or alkaline cleaners. Alcohol-based cleaning products are suitable. Maximum recommended cleaning temperature is 50°C. This product may be cleaned and reused by other persons.

ACCESSORIES AVAILABLE

- Wall strap kit - AWK
- Other products available include platforms.



Protected by UK patent applications and (European) Registered Community Designs. We undertake a policy of continual review and research, which may lead to changes in product appearance and specification. Date of Issue: October 2026.

Enable Access, Unit R, Cherrycourt Way, Leighton Buzzard, LU7 4UH, UK.

EC REP Enable Access Global Ltd, 27 Phibsboro Place, Dublin 7, D07 V20, Republic of Ireland.

For technical support or details on your nearest stockist please call 020 8275 0375.

These instructions are also available in an accessible PDF format suitable for blind and visually impaired persons.




Code	Length	Width	Lip Width	Capacity	Weight	Stored Dimensions L x W x D
AXM06	60cm (2ft)	79cm (31ins)	69.6cm (27.4ins)	525kg (82st)	5.91kg (13lbs)	60 x 43 x 10cm (24 x 17 x 4ins)
AXM09	90cm (3ft)	79cm (31ins)	69.6cm (27.4ins)	525kg (82st)	7.35kg (16.2lbs)	90 x 43 x 10cm (35 x 17 x 4ins)
AXM12	120cm (4ft)	79cm (31ins)	69.6cm (27.4ins)	475kg (74st)	10.63kg (23.4lbs)	120 x 43 x 10cm (47 x 17 x 4ins)
AXM15	150cm (5ft)	79cm (31ins)	69.6cm (27.4ins)	425kg (66st)	11.45kg (25.24lbs)	150 x 43 x 10cm (59 x 17 x 4ins)
AXM18	180cm (6ft)	79cm (31ins)	69.6cm (27.4ins)	425kg (66st)	15.46kg (34.08lbs)	180 x 43 x 10cm (71 x 17 x 4ins)
AXM21	210cm (7ft)	79cm (31ins)	69.6cm (27.4ins)	375kg (59st)	17.85kg (39.35lbs)	210 x 43 x 10cm (83 x 17 x 4ins)
AXM24	240cm (8ft)	79cm (31ins)	69.6cm (27.4ins)	325kg (51st)	20.24kg (44.62lbs)	240 x 43 x 10cm (95 x 17 x 4ins)

RAMP LENGTH



a. Measure the step height/level difference.
(Total height of steps or height to resting point on vehicle).

b. Use the **Ramp Length Guide**

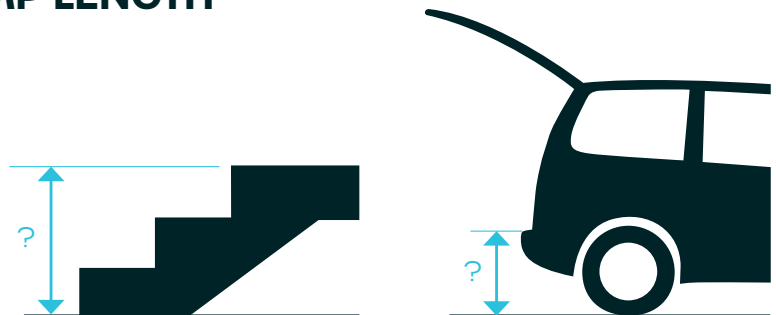

	Height (up to)	Length		
		1:6	1:8	1:12
Ramp Length Guide: Height x Gradient = Length Required	10cm (4ins)	60cm (2ft)	80cm (2ft 8ins)	120cm (4ft)
	15cm (6ins)	90cm (3ft)	120cm (4ft)	180cm (6ft)
	20cm (8ins)	120cm (4ft)	160cm (5ft 4ins)	240cm (8ft)
	25cm (10ins)	150cm (5ft)	200cm (6ft 8ins)	300cm (10ft)
	30cm (1ft)	180cm (6ft)	240cm (8ft)	360cm (12ft)
	45cm (1ft 6ins)	270cm (9ft)	360cm (12ft)	540cm (18ft)
	60cm (2ft)	360cm (12ft)	480cm (16ft)	720cm (24ft)

GRADIENT SLOPE

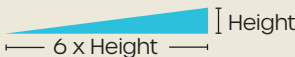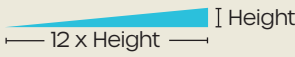1:6 (10°)*		Acceptable in some cases for short threshold ramps.
1:8 (7°)*		A useful balance between gradient and space suitable for portable ramps in some assisted applications and events, etc.
1:12 (5°)		Often required by social services and local authorities. Minimum required by Building Regulations Document M for installed ramps.
1:15 (4°)		Sometimes required by social services and local authorities for independent access but typically for installed public access ramps.

*Assess safety/risk