Enable Access

USER INSTRUCTIONS

Ultility Ramp

For products: UR9, UR12, UR15, UR18, UR21, UR24, URW9, URW12, URW15, URW18, URW21, URW24

BEFORE USE

- Remove all packaging and ensure the ramp has not been damaged in transit. Retain packaging.
- Carefully read the instructions and warning texts before using the ramps
- Keep instructions for reference
- On the ramp you will find a label that shows the maximum permitted load for the ramp.

Scan for more product information

TO DEPLOY

1. Carry/move the ramp using the handle and wheels (where fitted – see table page 3).
2. Ensure the ground/surface at both upper and lower levels is firm and stable.
3. Check ramp is long enough to cover the step/rise. The ramp is designed for use between 1:4 and 1:20 gradient (dependent on safe application and use requirements).
4. Place the ramp in position with the lip on the upper level, central to the chosen route.
5. Ensure the ramps do not become dislodged in use. Check ramp is firmly located and stable each time it is used.



MAINTENANCE

If the ramp becomes damaged or slippery they should not be used. The ramps may be cleaned with normal cleaning products. Maximum recommended cleaning temperature: 50°C.

RECYCLING

If you no longer require this product, please liaise with a metals recycling contractor. If this product was loaned to you by a local authority, please contact them to discuss its removal.

WARRANTY

- The ramp carries a comprehensive 2 year warranty for normal use, conditions apply.
- Use other than published/intended will void the warranty.
- These ramps are not designed for permanent outside use and must not be fixed down.
- Any modification/repair other than by manufacturer will void the warranty.
- Ramps will not be repaired/replaced if the serial number plaque has been removed/replaced

CE MD

Protected by UK patent applications and (European) Registered Community Designs. We undertake a policy of continual review and research, which may lead to changes in product appearance and specification. Date of Issue: July 2025.

Enable Access, Unit R, Cherrycourt Way, Leighton Buzzard, LU7 4UH, UK.

EC REP Enable Access Global Ltd, 27 Phibsboro Place, Dublin 7, D07 V20, Republic of Ireland.

For technical support or details on your nearest stockist please call 020 8275 0375.

MADE IN BRITAIN

Enable Access

USER INSTRUCTIONS

Utility Ramp

For products: UR9, UR12, UR15, UR18, UR21, UR24, URW9, URW12, URW15, URW18, URW21, URW24

UTILITYRAMP HANDLE KIT AND UTILITYRAMP WHEEL KIT FITTING INSTRUCTIONS:

- All UR and URW ramps use the Utilityramp Handle Kit (URHK)
- UR and URW ramps 150cm and longer also use the Utilityramp Wheel Kit (URWK)
- Casters are not fitted to UR9, UR12, URW9, URW12


Handle kit (URHK) contents:

- HD Ramp Peg (UR) x2
- M5 x 16 Coach Bolt x2
- Handle x1
- M5 Flange Nuts x2
- M8 x 55 Coach Bolt x2
- M8 Flange Nuts x2
- 1/4" Tommy Bar
- 7/8" Tube Spanner
- 13mm Tube Spanner

To fit handle:

1. Locate 2 x pre-drilled holes on the side of the ramp
2. Fit handle using M5 x 16mm coach bolt x2 and M5 Flange nut x2 (fixing tools supplied)



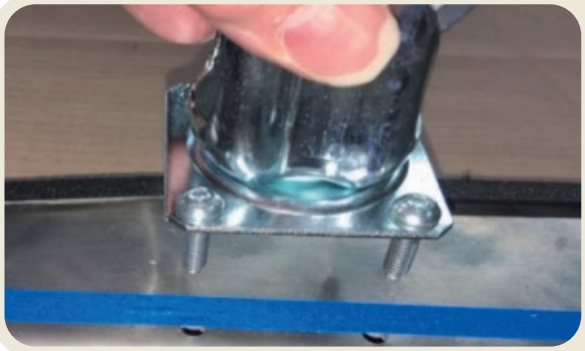

Wheel kit (URWK) contents:

- M6 x 20mm Coach Bolt x8
- M6 Flange nut x8
- 60mm Wheel Fixed Plate Castor x2
- 10mm Tube Spanner

To fit castors:

1. Locate the castors to two sets of 4x holes on the side of the ramp
2. Fix with M5 coach bolts x8 and M5 flange nut x8 (fixing tools supplied)

To fit ramp pegs (ARLS):

1. Locate the ramp pegs to the underside of the ramp at the threshold end on the 2nd set of holes from the end of the ramp.
2. Fix with M8 x 50 coach bolts x2 and M8 flange nyloc nuts x 2 using T bar and socket provided in kit (a screwdriver will also be required)



MAINTENANCE

If the ramp becomes damaged or slippery they should not be used. The ramps may be cleaned with normal cleaning products. Maximum recommended cleaning temperature: 50°C.

RECYCLING

If you no longer require this product, please liaise with a metals recycling contractor. If this product was loaned to you by a local authority, please contact them to discuss its removal.

WARRANTY

- The ramp carries a comprehensive 2 year warranty for normal use, conditions apply.
- Use other than published/intended will void the warranty.
- Any modification/repair other than by manufacturer will void the warranty.
- Ramps will not be repaired/replaced if the serial number plaque has been removed/replaced.

CE MD

Protected by UK patent applications and (European) Registered Community Designs. We undertake a policy of continual review and research, which may lead to changes in product appearance and specification. Date of Issue: July 2025.

Enable Access, Unit R, Cherrycourt Way, Leighton Buzzard, LU7 4UH, UK.

EC REP Enable Access Global Ltd, 27 Phibsboro Place, Dublin 7, D07 V20, Republic of Ireland.

For technical support or details on your nearest stockist please call 020 8275 0375.

MADE IN BRITAIN

Enable Access

USER INSTRUCTIONS

Ultility Ramp


Code	Length	Width	Capacity	Weight
UR9	90cm (3ft)	76cm (30ins)	400kg (63st)	9kg (20lbs)
UR12	120cm (4ft)	76cm (30ins)	400kg (63st)	12.5kg (28lbs)
UR15	150cm (5ft)	76cm (30ins)	400kg (63st)	15.5kg (34lbs)
UR18	180cm (6ft)	76cm (30ins)	400kg (63st)	19kg (42lbs)
UR21	210cm (7ft)	76cm (30ins)	350kg (55st)	22.5kg (50lbs)
UR24	240cm (8ft)	76cm (30ins)	300kg (47st)	25.5kg (56lbs)
URW9	90cm (3ft)	90cm (35ins)	400kg (63st)	10.5kg (23lbs)
URW12	120cm (4ft)	90cm (35ins)	400kg (63st)	14.5kg (32lbs)
URW15	150cm (5ft)	90cm (35ins)	400kg (63st)	18.5kg (41lbs)
URW18	180cm (6ft)	90cm (35ins)	400kg (63st)	22kg (49lbs)
URW21	210cm (7ft)	90cm (35ins)	350kg (55st)	26kg (57lbs)
URW24	240cm (8ft)	90cm (35ins)	300kg (47st)	30kg (66lbs)

RAMP LENGTH

a. Measure the step height/level difference.
(Total height of steps or height to resting point on vehicle).

b. Use the **Ramp Length Guide**



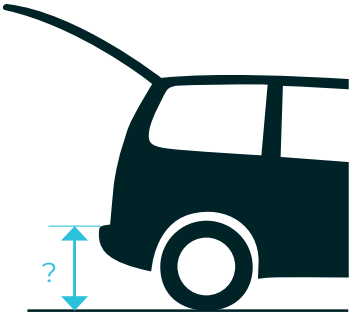

	Height (up to)	Length		
		1:6	1:8	1:12
Ramp Length Guide: Height x Gradient = Length Required	10cm (4ins)	60cm (2ft)	80cm (2ft 8ins)	120cm (4ft)
	15cm (6ins)	90cm (3ft)	120cm (4ft)	180cm (6ft)
	20cm (8ins)	120cm (4ft)	160cm (5ft 4ins)	240cm (8ft)
	25cm (10ins)	150cm (5ft)	200cm (6ft 8ins)	300cm (10ft)
	30cm (1ft)	180cm (6ft)	240cm (8ft)	360cm (12ft)
	45cm (1ft 6ins)	270cm (9ft)	360cm (12ft)	540cm (18ft)
	60cm (2ft)	360cm (12ft)	480cm (16ft)	720cm (24ft)

GRADIENT SLOPE

1:6 (10˚)*	Height / 6 x Height	Acceptable in some cases for short threshold ramps.
1:8 (7˚)*	Height / 8 x Height	A useful balance between gradient and space suitable for portable ramps in some assisted applications and events, etc.
1:12 (5˚)	Height / 12 x Height	Often required by social services and local authorities. Minimum required by Building Regulations Document M for installed ramps.
1:15 (4˚)	Height / 15 x Height	Sometimes required by social services and local authorities for independent access but typically for installed public access ramps.

*Assess safety/risk

CE MD

Protected by UK patent applications and (European) Registered Community Designs. We undertake a policy of continual review and research, which may lead to changes in product appearance and specification. Date of Issue: July 2025.

Enable Access, Unit R, Cherrycourt Way, Leighton Buzzard, LU7 4UH, UK.
EC REP Enable Access Global Ltd, 27 Phibsboro Place, Dublin 7, D07 V20, Republic of Ireland.

For technical support or details on your nearest stockist please call 020 8275 0375.

MADE IN BRITAIN